



Transfer spike and Draw-off tube solutions

**Sterimatic Worldwide Ltd
Griffin Mill, GL5 2AZ, UK.**

**www.sterimatic.com
info@sterimatic.com**


Tel: 0044 1453 884944



Transfer Spike solutions from Sterimatic

We have several devices that can be used to safely mix fluids with fluids or fluids with powders. These devices can be supplied individually packed and sterile or in bulk.

Standard Transfer Spike Code T22

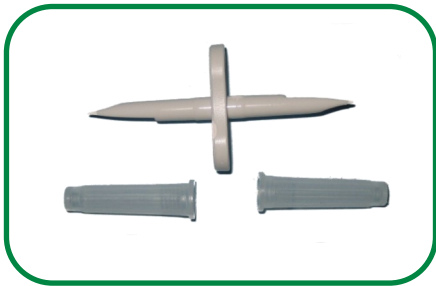


This plastic transfer spike can be used for safe mixing.

It is supplied with 2 end caps to protect the spike in transit.

Each device is individually packed and sterilised by EO.

Sterimatic Transfer Spike Code T01

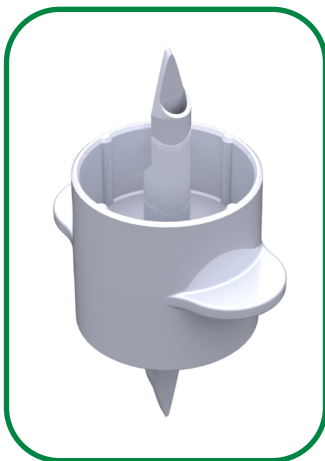

We have developed a new transfer spike;

- Fast transfer of fluid
- Strong spike
- Support cups to hold the spike onto the bottle.
- Empties the vials
- Suitable for 20mm diameter necks.
- Customisable colour

Can be supplied Sterile or non-sterile



Spike + draw-off tube solutions from Sterimatic

We have several devices that can be used to withdraw fluids from bottles or collapsible bags. These kits can be customised to meet your requirements with bespoke components, tube size and lengths, printed inserts, ancillary items etc.....
Below is a selection of kits. please contact us if you require something different.

Vented Spike with Draw-off tube



One solution for all multi-dose tube-fed bottles.


The air filter reduces bottle contamination and the closure cap enables it to be used with standard bottles or collapsible packs. The spike fits securely in the bottle stopper, this enables the same device to be used for different bottle neck sizes.

Options:

1. **SYX012BS:** Vented spike + 1m draw-off tube (sterile)
2. **SYX015STE:** Vented spike + 1.5m draw-off tube (sterile)

Packs can be supplied sterile or non-sterile and include an insert.

Non-Vented Spike with Draw-off tube

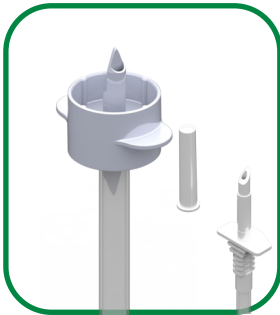


For use with collapsible bags only.

The non-vented spike allows fluid to withdraw from the bottle and as it does so the bag collapses.
2 options. Universal spike for all bottle necks or we have a new spike with cap for 20mm bottle necks.

1. **SYX007M:** Non-Vented universal spike + 1m draw-off tube (sterile).
2. **SYX008M:** Non-Vented spike (20mm) + 1m draw-off tube (sterile)

Packs can be supplied sterile or non-sterile and include an insert.

Spikes with 2 (or 4) Draw-off tubes

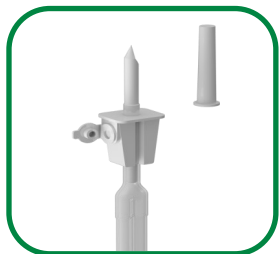

Where 2 operators will use the same bottle of vaccine

Options:

- Vented spike + 2 x 1.6m draw-off tubes (sterile)
- Non-Vented spike + 2 x 1.6m draw-off tubes (sterile)

Packs can be supplied sterile or non-sterile and include an MSD insert.

Bottle Clip. Code: CLP001



The spike holds securely inside the rubber stopper, however if you require extra security, we have a bottle clip available for 20 & 27mm diameter necked bottles.

Simply insert the spike and then fit the clip around the spike and bottle neck.




IV Kit solutions from Sterimatic

For small doses where it is important to see the flow of vaccine a drip chamber can be added. We can also include a tube clamp to enable changing the spike between bottles without getting air-bubbles into the tubing.

These kits can be customised to meet your requirements with bespoke components, tube size and lengths, printed inserts, ancillary items etc.

Below is a selection of available kits however please contact us if you require something different.

Vented Spike & Drip chamber with Draw-off tube

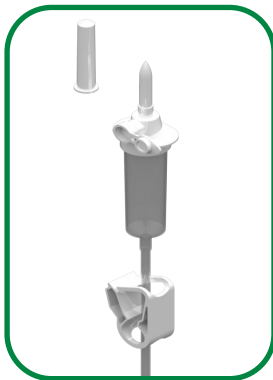


The kit can be supplied on its own or with a needle & syringe for mixing concentrates and diluents

Options:

- **SYA015:** Vented spike & drip chamber + 0.6m (or 0.8 or 1.0m) draw-off tube (4.5mm ID) + tube clamp
- **SYA0010I:** As above with 5ml disposable syringe and 18G x 1 1/2" disposable Needle.

Supplied Sterile and includes an insert.

Non-Vented Spike & Drip chamber with Draw-off tube

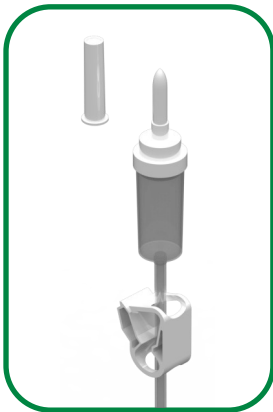

For use with collapsible bags only.

The non-vented spike allows fluid to pass from the bottle and as it does so the bag collapses.


The spike fits securely in the bottle stopper, this enables the same device to be used for different bottle neck sizes.

Options:

- **SYA095:** Vented spike & drip chamber + 0.6m (or 0.8 or 1.0m) draw-off tube (4.5mm ID) + tube clamp
- **SYA096M:** As above with 5ml disposable syringe and 18G x 1 1/2" Needle

Supplied Sterile and includes an insert.

Spike & Drip chamber with 2 Draw-off tubes

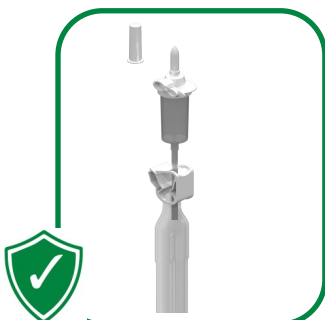

Where 2 operators will use the same bottle of vaccine

Same options as above. Either vented or non-vented spike & drip chamber + tubing.



Contact details:

Sterimatic Worldwide Ltd
Griffin Mill, GL5 2AZ, UK.

www.sterimatic.com
info@sterimatic.com

Tel: 0044 1453 884944



Certificate Number: 16127
ISO 9001



Sterimatic products are protected by patents, design registrations, trademarks and copyrights.
This brochure and its contents are copyrighted by Sterimatic Worldwide Ltd 2018 and all rights are reserved.

LEF041 04/2018